


Technical Information about Sp-ITP- α -S

Update: March 17, 2023ss

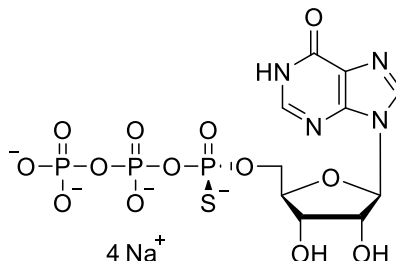

Abbreviation: Sp-ITP- α -S

Formula	CAS No.	Molecular Weight	UV	BIOLOG Cat. No.
C ₁₀ H ₁₅ N ₄ O ₁₃ P ₃ S (free acid)	[347887-30-3]	524.2 (free acid)	λ_{\max} 249 nm / ϵ 12800 / pH7	I 026

Name: Inosine- 5'- O- (1- thiotriphosphate), Sp- isomer (Sp-ITP- α -S), sodium salt

Description: Sp-ITP- α -S is a modification of inosine triphosphat (ITP), where one of the non-bridging oxygens in the S position of the α - phosphate is replaced by sulfur. The suffix "p" indicates that R/S nomenclature refers to phosphorus. The corresponding Rp-isomer is offered as well (Cat. No. I 027).

Properties:

- Increased metabolic stability compared to ITP,
- for research into ATP-dependent receptor proteins and for determination of their stereospecificity,
- accepted by RNA polymerase such as T7 RNA polymerase for incorporation of phosphorothioate into RNA.

Specification: Aqueous solution of the sodium salt (10 mM). Other salts of Sp-ITP- α -S are available upon request. Micromolar quantities are determined by UV at λ_{\max} . When opening the tube please make sure that no liquid is lost within the cap. A short spin-down in a bench centrifuge is recommended before use.

Purity: Typical analysis is better than 95% (HPLC / UV / 249 nm) at time of quality control and packing. However, actual purity depends on storage and transport conditions. The product is not sterile and has not been tested for endotoxins.

Stability and Storage: Sp-ITP- α -S is most stable when stored as aqueous solution in the freezer (-20° Celsius necessary, -70° recommended), however, at ambient temperature the compound slowly starts to decompose. Thus, in order to maintain its original high quality it is recommended to allow thawing only before using the product. If you will not use up the vial with one application, please aliquot the contents of the vial in order to avoid repeated freeze/thaw cycles for the rest. When making such aliquots be sure to operate quickly and to freeze the vial again as soon as possible.

Toxicity and Safety: Since nucleoside triphosphates have multiple tasks in every organism, it is very likely that triphosphate analogues will interfere with many cell regulation processes *in vivo*. However, due to the rather small quantities to work with, no health hazards have been reported. Nevertheless please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid skin contact or ingestion and allow only trained personnel to handle the product.

Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

Not for drug, household or related uses!

Selected Reference for Sp-ITP- α -S:

Heide, C.; Busch, S.; Feltens, R.; Hartmann, R. K., *RNA*, **7**, 553 - 564 (2001): "Distinct Modes of Mature and Precursor tRNA Binding to *Escherichia coli* RNase P RNA Revealed by NAIM Analyses"

Selected References for the Related Compound Sp-ATP- α -S (Cat. No. A 040):

- Hillaire-Buys, D.; Shahar, L.; Fischer, B.; Chulkin, A.; Linck, N.; Chapal, J.; Loubatieres-Mariani, M. M.; Petit, P., *Drug Dev. Res.*, **53**, 33 - 43 (2001): "Pharmacological Evaluation and Chemical Stability of 2-Benzylthioether-5'-O-(1-Thiotriphosphate)-Adenosine, a new Insulin Secretagogue Acting Through P2Y Receptors"
- Laubinger, W.; Reiser, G., *Biochem. Pharmacol.*, **55**, 687 - 695 (1998): "Differential Characterization of Binding Sites for Adenine and Uridine Nucleotides in Membranes from Rat Lung as Possible Tools for Studying P2 Receptors in Lung"
- Conigrave, A.D.; Lee, J.Y.; van der Weyden, L.; Jiang, L.; Ward, P.; Tasevski, V.; Luttrell, B.M.; Morris, M.B., *Br. J. Pharmacol.*, **124**, 1580 - 1585 (1998): "Pharmacological Profile of a Novel Cyclic AMP-linked P2 Receptor on Undifferentiated HL-60 Leukemia Cells"
- Xu, S.H.; Gaskin, F., *Biochim. Biophys. Acta*, **1383**, 111 - 122 (1998): "Probing the ATP Binding Site of Tubulin with Thiotriphosphate Analogues of ATP"
- Hechler, B.; Leon, C.; Vial, C.; Vigne, P.; Cazenave, J.P.; Gachet, C., *Blood*, **92**, 152 - 159 (1998): "The P2Y(1) Receptor is Necessary for Adenosine 5'-diphosphate-induced Platelet Aggregation"
- Vörtler, C.S.; Fedorova, O.; Persson, T.; Kutzke, U.; Eckstein, F., *RNA*, **4**, 1444 - 1454 (1998): "Determination of 2'-Hydroxyl and Phosphate Groups Important for Aminoacylation of Escherichia coli tRNA Asp: A Nucleotide Analogue Interference Study"
- King, B.F.; Wildmann, S.S.; Ziganshina, L.E.; Pintor, J.; Burnstock, G., *Br. J. Pharmacol.*, **121**, 1445 - 1453 (1997): "Effects of Extracellular pH on Agonism and Antagonism at a Recombinant P2X(2) Receptor"
- Mitchell, C.; Blaho, J.A.; McCormick, A.L.; Roizman, B., *J. Biol. Chem.*, **272**, 25394 - 25400 (1997): "The Nucleotidylylation of Herpes Simplex virus 1 Regulatory Protein Alpha-22 by Human Casein Kinase II"
- Cullis, P.M., Maxwell, A.; Weiner, D.P., *Biochemistry*, **36**, 6059 - 6068 (1997): "Exploiting Nucleotide Thiophosphates to Probe Mechanistic Aspects of Escherichia coli DNA Gyrase"
- Cusack, N.J.; Hourani, S.M.O., *Nucleosides Nucleotides*, **10**, 1019 - 1028 (1991): "Design, Syntheses and Pharmacology of ATP Analogues Selective for Subtypes of P2-Purinoceptors"
- Lazewska, D.; Guranowski, A., *Nucleic Acids Res.*, **18**, 6083 - 6088 (1990): "P- α -chiral Phosphorothioate Analogues of bis (5'-adenosyl)tetraphosphate (Ap4A; their Enzymatic Synthesis and Degradation"
- Eckstein, F.; Gish, G., *TIBS*, **14**, 97 - 100 (1989): "Phosphorothioates in Molecular Biology"
- Griffiths, A.D.; Potter, B.V.L.; Eperon, I.C., *Nucl. Acids Res.*, **15**, 4145 - 4162 (1987): "Stereospecificity of Nucleases Towards Phosphorothioate-substituted RNA: Stereochemistry of Transcription by T7 RNA Polymerase"
- Eckstein, F., *Ann. Rev. Biochem.*, **54**, 367 - 402 (1985): "Nucleoside Phosphorothioates"
- Shigalowa, T.; Lehmann, U.; Krevet, M.; Strotmann, H., *Biochim. Biophys. Acta*, **809**, 57 - 65 (1985): "Transient Stimulation of Light-triggered Hydrolysis by Preillumination of Chloroplasts in the Presence of ATP"
- Cusack, N.J.; Pearson, J.D.; Gordon, J.L., *Biochem. J.*, **214**, 975 - 981 (1983): "Stereoselectivity of Ectonucleotidases on Vascular Endothelial Cells"
- Eckstein, F., *Angew. Chemie*, **95**, 431 - 451 (1983): "Phosphorothioatanaloga von Nucleotiden - Werkzeuge zur Untersuchung biochemischer Prozesse"
- Marlier, J.F.; Benkovic, S.J., *Biochemistry*, **21**, 2349 - 2356 (1982): "On the Mechanism of de Novo Polymerization by Form I Polynucleotide Phosphorylase of *Micrococcus luteus*"
- Cusack, N.J.; Hourani, S.M.O., *Br. J. Pharmacol.*, **76**, 221 - 227 (1982): "Adenosine 5'-Diphosphate Antagonists and Human Platelets: No Evidence that Aggregation and Inhibition of Stimulated Adenylate Cyclase are Mediated by Different Receptors"
- Lowe, G.; Tansley, G.; Cullis, P.M., *J. Chem. Soc. Chem. Commun.*, **595** (1982): "Synthesis of Adenosine- 5' [(R)alpha-¹⁷O] Triphosphate"
- Marlier, J.F.; Bryant, F.R.; Benkovic, S.J., *Biochemistry*, **20**, 2212 - 2219 (1981): "Stereochemical and Kinetic Investigation of 32P-Labeled Inorganic Phosphate Exchange Reaction Catalyzed by Primer-independent and Primer-dependent Polynucleotide Phosphorylase from *Micrococcus luteus*"
- Eckstein, F.; Romaniuk, P.J.; Heidemann, W.; Storm, D.R., *J. Biol. Chem.*, **256**, 9118 - 9129 (1981): "Stereochemistry of the Mammalian Adenylate Cyclase Reaction"
- Strotmann, H.; Bickel-Sandkötter, S.; Edelman, K.; Eckstein, F.; Schlimme, E.; Boos, K.S.; Lüstorf, J., *Biochim. Biophys. Acta*, **545**, 122 - 130 (1979): "Thiophosphate Analogs of ADP and ATP as Substrates in Partial Reactions of Energy Conversion in Chloroplasts"

Ngoc, L.D.; Jebeleanu, G.; Barzu, O., *FEBS Lett.*, **97**, 65 - 68 (1979): "Interaction of Sulfur-containing ATP Analogs With Rabbit Muscle Phosphofructokinase"

Yee, D., Armstrong, V.W.; Eckstein, F., *Biochem.*, **18**, 4116 - 4120 (1979): "Mechanistic Studies on Deoxyribonucleic Acid Dependent Ribonucleic Acid Polymerase from Escherichia coli Using Phosphorothioate Analogues. 1. Initiation and Pyrophosphate Exchange Reactions"

Burgers, P.M.J.; Eckstein, F., *Proc. Natl. Acad. Sci. USA*, **75**, 4798 - 4800 (1978): "Absolute Configuration of the Diastereomers of Adenosine 5'-O-(1-thiotriphosphate): Consequences for the Stereochemistry of Polymerization by DNA-dependent RNA Polymerase from Escherichia coli"