


## Technical Information about 2',3'-cAMP

Update September 29, 2022 ss

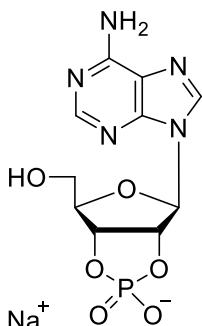

**Abbreviation:** 2',3'-cAMP

Formula	CAS No.	Molecular Weight	UV	BIOLOG Cat. No.
C <sub>10</sub> H <sub>11</sub> N <sub>5</sub> O <sub>6</sub> P·Na	[37063-35-7]	351.2	λ <sub>max</sub> 259 nm / ε 15000 / pH 7	A 307

**Name:** Adenosine- 2', 3'- cyclic monophosphate, sodium salt

**Description:** 2',3'-cAMP is a cyclic phosphate ester of adenosine in which both the 2'- and the 3'-hydroxy groups are esterified by phosphoric acid.

**Properties:** 2',3'-cAMP is a positional isomer of the natural second messenger 3',5'-cAMP (BIOLOG Cat. No. A 001) and was found to be produced by organ systems such as rat and mouse kidney and mouse brain from RNA degradation. 2',3'-cAMP is metabolized to 2'-AMP and 3'-AMP, and these monophosphates are subsequently converted to adenosine. It was also recently reported to promote the formation of stress granules in *Arabidopsis thaliana*. Several pathways regulated by 2',3'-cNMPs, including flagellar motility and biofilm formation in *E. coli* were identified by using the prodrug / precursor 5'-O-MB-2',3'-cAMP (BIOLOG Cat. No. M 176).

**Specification:** Lyophilized or crystallized sodium salt. Other salt forms are available upon request. Equal concentrations of 2',3'-cAMP can appear very different in volume due to sensitivity of the lyophilized form to humidity. The compound can even contract to small volume droplets. Normally the product is located in the conical bottom of the tube. Micromolar quantities are determined by UV at λ<sub>max</sub>.

**Purity:** Typical analysis is better than 97 % (HPLC / UV / 259 nm). The product is not sterile and has not been tested for endotoxins.

**Solubility:** 2',3'-cAMP is soluble in water (≥ 350 mM, limits have not been determined). Please rinse tube walls carefully and preferably use ultrasonic or vortex to achieve total and uniform mixing. When opening the tube please make sure that no substance is lost within the cap.

**Stability and Storage:** 2',3'-cAMP has sufficient stability at room temperature and does not need special care during handling or shipment. Nevertheless, we recommend that the compound should be stored in the freezer, for longer storage periods preferably in freeze-dried form.

**Toxicity and Safety:** Please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid skin contact or ingestion and allow only trained personnel to handle the product. Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

**Not for drug, household or related uses!**

### Selected References for 2',3'-cAMP:

Duggal, Y.; Kurasz, J. E.; Fontaine, B. M.; Marotta, N. J.; Chauhan, S. S.; Karls, A. C.; Weinert, E. E., *J. Bacteriol.*, **204**, e0020821 (2022): "Cellular Effects of 2',3'-Cyclic Nucleotide Monophosphates in Gram-Negative Bacteria"

Ludwig, N.; Yerneni, S.S.; Menshikova, E.V.; Gillespie, D.G.; Jackson, E.K.; Whiteside, T.L., *Scientific Reports*, **10**, 6948 (2020): "Simultaneous Inhibition of Glycolysis and Oxidative Phosphorylation Triggers a Multi-Fold Increase in Secretion of Exosomes: Possible Role of 2'3'-cAMP"

Kosmacz, M.; Luzarowski, M.; Kerber, O.; Leniak, E.; Gutiérrez-Beltrán, E.; Moreno, J.C.; Gorka, M.; Szlachetko, J.; Veyel, D.; Graf, A.; Skirycz, A., *Plant Physiol.*, **177**, 411 - 421 (2018): "Interaction of 2',3'-cAMP with Rbp47b Plays a Role in Stress Granule Formation"

Forman, M.B.; Gillespie, D.G.; Cheng, D.; Jackson, E.K., *Dig. Dis. Sci.*, **59**, 2118 - 2125 (2014): "A Novel Adenosine Precursor 2',3'-Cyclic Adenosine Monophosphate Inhibits Formation of Post-Surgical Adhesions"