


Technical Information about pA-3'-p / pAp

Update: January 26, 2024 ss

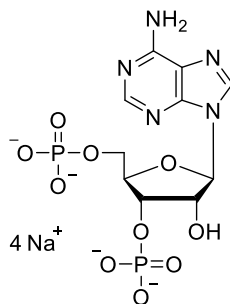

Abbreviation:

pA-3'-p

Formula	CAS No.	Molecular Weight	UV	BIOLOG Cat. No.
C ₁₀ H ₁₅ N ₅ O ₁₀ P ₂ (free acid)	[75431-54-8]	427.2 (free acid)	λ _{max} 259 nm / ε 15000 / pH 7	A 444

Name: Adenosine- 3', 5'- O- bisphosphate, sodium salt

Syn.: PAP, 3'-PAP / 3'- phosphoadenosine 5'- phosphate / adenosine- 3', 5'- diphosphate

Description: In pA-3'-p phosphate groups are attached to both the ribose 3' position and the ribose 5' position of adenosine.

Properties: pA-3'-p is an ubiquitous nucleotide present in almost all organisms. It is a product of 3'-phosphoadenosine 5'-phosphosulfate (PAPS), a cofactor of sulfotransferases (ST / SULTs), and can be used to study the kinetic properties and structure of sulfotransferases. pA-3'-p can also act as an inhibitor of exoribonucleases, such as Xrn1, Rat1p and DXO (Dichtel et al. 1997, Yun et al. 2018), and of poly(ADP-ribose) polymerase 1 (PARP-1) (Toledano et al. 2012). Furthermore, it inhibits ADP-induced platelet shape change and aggregation but it does not affect ADP-induced inhibition of adenylyl cyclase (Hechler et al. 1998).

Specification: Crystallized or lyophilized sodium salt. For other salt forms or analogues of pA-3'-p please inquire. Please keep in mind that equal concentrations of the compound may look different in volume due to high sensitivity of the lyophilized form to humidity. The compound can even contract to small volume droplets. Normally the product is located in the conical bottom of the tube. Micromolar quantities are determined by UV at λ_{max}.

Purity: Typical analysis is better than 95% (HPLC / UV / 259 nm). The product is not sterile and has not been tested for endotoxins.

Solubility: pA-3'-p is soluble in water (≥ 38 mM, limits have not been determined). Please rinse tube walls carefully and preferably use ultrasonic or vortex to achieve total and uniform mixing. When opening the tube please make sure that no substance is lost within the cap.

Stability and Storage: pA-3'-p is chemically rather stable and does not need special care during handling or shipment. Nevertheless, we recommend that the compound should be stored in the freezer, for longer storage periods preferably in freeze-dried form.

Toxicity and Safety: Please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid skin contact or ingestion and allow only trained personnel to handle the product.

Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

Not for drug, household or related uses!

Selected References for pA-3'-p:

Yun, J. S.; Yoon, J. H.; Choi, Y. J.; Son, Y. J.; Kim, S.; Tong, L.; Chang, J. H., *Biochem. Biophys. Res. Commun.*, **504**, 89 - 95 (2018): "Molecular Mechanism for the Inhibition of DXO by adenosine 3',5'-bisphosphate"

Tibbs, Z. E.; Rohn-Glowacki, K. J.; Crittenden, F.; Guidry, A. L.; Falany, C. N., *Drug Metab. Pharmacokinet.*, **30**, 3 - 20, (2015): "Structural Plasticity in the Human Cytosolic Sulfotransferase Dimer and its Role in Substrate Selectivity and Catalysis"

Toledano, E.; Ogryzko, V.; Danchin, A.; Ladant, D.; Mechold, U., *Biochem. J.*, **443**, 485 - 490 (2012): "3'-5' Phosphoadenosine Phosphate is an Inhibitor of PARP-1 and a Potential Mediator of the Lithium-Dependent Inhibition of PARP-1 in vivo"

Estavillo, G.M.; Crisp, P.A.; Pornsiriwong, W.; Wirtz, M.; Collinge, D.; Carrie, C.; Giraud, E.; Whelan, J.; David, P.; Javot, H.; Brearley, C.; Hell, R.; Marin, E.; Pogson, B.J., *Plant Cell*, **23**, 3992 - 4012 (2011): "Evidence for a SAL1-PAP Chloroplast Retrograde Pathway That Functions in Drought and High Light Signaling in *Arabidopsis*"

Dombrowski, L.; Dong, A.; Bochkarev, A.; Plotnikov, A. N., *Proteins*, **64**, 1091 - 1094 (2006): "Crystal Structures of Human Sulfotransferases SULT1B1 and SULT1C1 Complexed with the Cofactor Product adenosine-3'-5'-diphosphate (PAP)"

Sheng, J. J.; Duffel, M. W., *Protein Expr. Purif.*, **21**, 235 - 242 (2001): "Bacterial Expression, Purification, and Characterization of Rat Hydroxysteroid Sulfotransferase STa"

Hechler, B.; Léon, C.; Vial, C.; Vigne, P.; Frelin, C.; Cazenave, J. P.; Gachet, C., *Blood*, **92**, 152 - 159 (1998): "The P2Y₁ Receptor is Necessary for Adenosine 5'-diphosphate-Induced Platelet Aggregation"

Dichtl, B.; Stevens, A.; Tollervey, D., *EMBO J.*, **16**, 7184 - 7195 (1997): "Lithium Toxicity in Yeast is Due to the Inhibition of RNA Processing Enzymes"

Boyer, J. L.; Romero-Avila, T.; Schachter, J. B.; Harden, T.K., *Mol. Pharmacol.*, **50**, 1323 -1329 (1996): "Identification of Competitive Antagonists of the P2Y₁ Receptor"