


Technical Information about 3'-Cyclic ADP Ribose

Update: July 19, 2023 ss

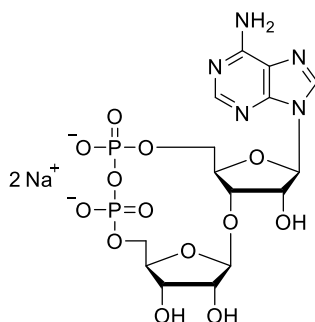

Abbreviation:

3'-cADPR

| Formula | CAS No. | Molecular Weight | UV | BIOLOG Cat. No. |
|--|-----------|----------------------|--|-----------------|
| C ₁₅ H ₂₁ N ₅ O ₁₃ P ₂ (free acid) | [pending] | 541.3 (free acid) | λ _{max} 259 nm / ε 15000 / pH 7 | C 404 |

Name: (1'' - 3')- Cyclic adenosine diphosphate ribose (3'-cADPR / v2-cADPR / 1''-3'-gcADPR), sodium salt

Description: 3'-cADPR is an analogue of ADP ribose in which position 3' of the southern ribose is connected to position 1'' of the second ribose unit yielding a so-called "glycocyclic" macrocyclic ring system consisting of pyrophosphate and riboses.

Properties: 3'-cADPR is a regioisomer of the canonical second messenger cADPR (Cat. No. C 005) and of 2'-cADPR (Cat. No. C 406). 3'-cADPR was identified as novel signaling molecule involved in the innate immunity in plants and bacteria (Leavitt et al. 2022, Manik et al. 2022).

Specification: Crystallized or lyophilized sodium salt. The free acid or other salt forms are available upon request. Please keep in mind that equal concentrations of the compound may look different in volume due to sensitivity of the lyophilized form to humidity. The compound can even contract to small volume droplets. Normally the product is located in the conical bottom of the tube. Micromolar quantities are determined by UV at λ_{max}.

Purity: Typical analysis is better than 98% (HPLC / UV / 259 nm). The product is not sterile and has not been tested for endotoxins.

Solubility: 3'-cADPR has excellent solubility in water and aqueous buffers (≥ 9.5 mM, limits have not been determined). Please rinse tube walls carefully and preferably use ultrasonic or vortex to achieve total and uniform mixing. When opening the tube please make sure that no substance is lost within the cap.

Stability and Storage: 3'-cADPR has sufficient stability at room temperature and does not need special care during handling or shipment. Nevertheless, we recommend that the compound should be stored in the freezer, for longer storage periods preferably in freeze-dried form.

Toxicity and Safety: Since 3'-cADPR seems to have tasks in several organisms, it is not unlikely that it will interfere with many cell regulation processes *in vivo*. However, due to the rather small quantities to work with, no health hazards have been reported. Nevertheless please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid skin contact or ingestion and allow only trained personnel to handle the product.

Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

Not for drug, household or related uses!

Selected References for 3'-cADPR:

Bayless, A.M.; Chen, S.; Ogden, S.C.; Xu, X.; Sidda, J.D.; Manik, M.K.; Li, S.; Kobe, B.; Ve, T.; Song, L.; Grant, M.; Wan, L.; Nishimura, M.T., *Sci. Adv.*, **9**, eade8487 (2023): " Plant and Prokaryotic TIR Domains Generate Distinct Cyclic ADPR NADase Products"

Hulin, M.T.; Hill, L.; Jones, J.D.G.; Ma, W., *Proc. Natl. Acad. Sci. U. S. A.*, **120**, e2217114120 (2023): "Pangenomic Analysis Reveals Plant NAD⁺ Manipulation as an Important Virulence Activity of Bacterial Pathogen Effectors"

Yirmiya, E.; Leavitt, A.; Lu, A.; Avraham, C.; Osterman, I.; Garb, J.; Antine, S.P.; Mooney, S.E.; Hobbs, S.J.; Kranzusch, P.J.; Amitai, G.; Sorek, R., *bioRxiv*, <https://doi.org/10.1101/2023.05.01.538930> (2023): "Phages Overcome Bacterial Immunity via Diverse Anti-Defense Proteins"

Eastman, S.; Bayless, A.; Guo, M., *Mol. Plant Microbe Interact.*, **35**, 964 - 976 (2022): "The Nucleotide Revolution: Immunity at the Intersection of Toll/Interleukin-1 Receptor Domains, Nucleotides, and Ca²⁺"

Eastman, S.; Smith, T.; Zaydman, M.A.; Kim, P.; Martinez, S.; Damaraju, N.; DiAntonio, A.; Milbrandt, J.; Clemente, T.E.; Alfano, J.R.; Guo, M., *New Phytol.*, **233**, 890 - 904 (2022): "A Phyto bacterial TIR Domain Effector Manipulates NAD⁺ to Promote Virulence"

Leavitt, A.; Yirmiya, E.; Amitai, G.; Lu, A.; Garb, J.; Herbst, E.; Morehouse, B.R.; Hobbs, S.J.; Antine, S.P.; Sun, Z.J.; Kranzusch, P.J.; Sorek, R., *Nature*, **611**, 326 - 331 (2022): "Viruses Inhibit TIR gcADPR Signalling to Overcome Bacterial Defence"

Manik, M.K.; Shi, Y.; Li, S.; Zaydman, M.A.; Damaraju, N.; Eastman, S.; Smith, T.G.; Gu, W.; Masic, V.; Mosaiab, T.; Weagley, J.S.; Hancock, S.J.; Vasquez, E.; Hartley-Tassell, L.; Kargios, N.; Maruta, N.; Lim, B.Y.J.; Burdett, H.; Landsberg, M.J.; Schembri, M.A.; Prokes, I.; Song, L.; Grant, M.; DiAntonio, A.; Nanson, J.D.; Guo, M.; Milbrandt, J.; Ve, T.; Kobe, B., *Science*, **377**, eadc8969 (2022): "Cyclic ADP Ribose Isomers: Production, Chemical Structures, and Immune Signaling"

Weagley, J.S.; Zaydman, M.; Venkatesh, S.; Sasaki, Y.; Damaraju, N.; Yenkin, A.; Buchser, W.; Rodionov, D.A.; Osterman, A.; Ahmed, T.; Barratt, M.J.; DiAntonio, A.; Milbrandt, J.; Gordon, J.I., *Cell Rep.*, **39**, 110738 (2022): "Products of Gut Microbial Toll/Interleukin-1 Receptor Domain NADase Activities in Gnotobiotic Mice and Bangladeshi Children with Malnutrition"

Ofir, G.; Herbst, E.; Baroz, M.; Cohen, D.; Millman, A.; Doron, S.; Tal, N.; Malheiro, D.B.A.; Malitsky, S.; Amitai, G.; Sorek, R., *Nature*, **600**, 116 - 120 (2021): "Antiviral Activity of Bacterial TIR Domains via Immune Signalling Molecules"

Wan, L.; Essuman, K.; Anderson, R.G.; Sasaki, Y.; Monteiro, F.; Chung, E.H.; Osborne Nishimura, E.; DiAntonio, A.; Milbrandt, J.; Dangl, J.L.; Nishimura, M.T., *Science*, **365**, 799 - 803 (2019): "TIR Domains of Plant Immune Receptors are NAD⁺-Cleaving Enzymes that Promote Cell Death"