


Technical Information about 2-MeS-ITP

Update: March 17, 2023 ss

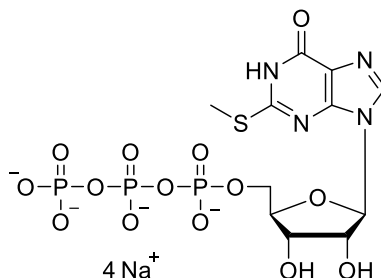

Abbreviation: 2-MeS-ITP

Formula	CAS No.	Molecular Weight	UV	BIOLOG Cat. No.
C ₁₁ H ₁₇ N ₄ O ₁₄ P ₃ S (free acid)	[pending]	554.3 (free acid)	λ _{max} 262 nm / ε 15400 / pH 7	M 114

Name: 2- Methylthioinosine- 5'- O- triphosphate (2-MeS-ITP), sodium salt

Description: 2-MeS-ITP is an analogue of 2-MeS-ATP (Cat. No. M 021) where the 6-amino group has been replaced by an oxygen atom.

Properties: 2-MeS-ITP is a potential metabolite of 2-MeS-ATP (Cat. No. M 021) and a potential agonist of P2Y purinoceptors. It is supposed to be an inhibitor of soluble guanylate cyclase.

Specification: Aqueous solution of the sodium salt (10 mM). Other salt forms of 2-MeS-ITP are available upon request. Micromolar quantities are determined by UV at λ_{max}. When opening the tube make sure that no liquid is lost within the cap. A short spin-down in a bench centrifuge is recommended before use.

Purity: Typical purity is better than 95% (HPLC / UV / 262 nm) at time of quality control and packing. However, actual purity depends on storage and transport conditions. The product is not sterile and has not been tested for endotoxins.

Stability and Storage: 2-MeS-ITP is most stable when stored as aqueous solution in the freezer (-20° Celsius necessary, -70° recommended), however, at ambient temperature the compound slowly starts to decompose. Thus, in order to maintain its original high quality it is recommended to allow thawing only before using the product. If you will not use up the vial with one application, please aliquot the contents of the vial in order to avoid repeated freeze/thaw cycles for the rest. When making such aliquots be sure to operate quickly and to freeze the vial again as soon as possible.

Toxicity and Safety: Since nucleoside triphosphates have multiple tasks in every organism, it is very likely that triphosphate analogues will interfere with many cell regulation processes *in vivo*. However, due to the rather small quantities to work with, no health hazards have been reported. Nevertheless please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid skin contact or ingestion and allow only trained personnel to handle the product.

Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

Not for drug, household or related uses!

Selected Reference for 2-MeS-ITP: 2-MeS-ITP is a new product and there are currently no corresponding references available.

Selected References for the Related Compound 2-MeS-ATP (Cat. No. M 021):

Swennen, E.L.; Bast, A.; Dagnelie, P.C., *Biochem. Biophys. Res. Commun.*, 348, 1194 - 1199 (2006): "Purinergic Receptors Involved in the Immunomodulatory Effects of ATP in Human Blood"

Ruiz-Stewart, I.; Kazerounian, S.; Pitari, G.M.; Schulz, S.; Waldman, S.A., *Eur. J. Biochem.*, **269**, 2186 - 2193 (2002): "Soluble Guanylate Cyclase is Allosterically Inhibited by Direct Interaction with 2-substituted Adenine Nucleotides"

Burnstock, G.; Cusack, N.J.; Hills, J.M.; MacKenzie, I.; Meghji, P., *Br. J. Pharmacol.*, **79**, 907 - 913 (1983): "Studies on the Stereoselectivity of the P2-Purinoceptor"

Cusack, N.J.; Hourani, S.M.O., *Br. J. Pharmacol.*, **75**, 397 - 400 (1982): "Specific but Noncompetitive Inhibition by 2-Alkylthio Analogues of Adenosine 5'-monophosphate and Adenosine 5'-triphosphate of Human Platelet Aggregation Induced by Adenosine 5'-diphosphate"