


## Technical Information about 2',2'',2'''-Tri-F-c-tridAMP

Fluorinated version of the bacterial second messenger c-triAMP

Update: February 19, 2021 AI

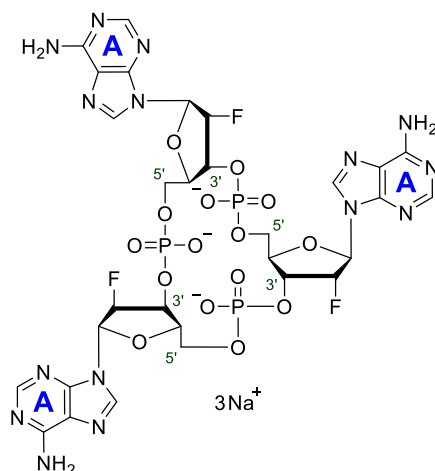

**Abbreviation:** 2',2'',2'''-Tri-F-c-tridAMP

Formula	CAS No.	Molecular Weight	UV	BIOLOG Cat. No.
C <sub>30</sub> H <sub>33</sub> F <sub>3</sub> N <sub>15</sub> O <sub>15</sub> P <sub>3</sub> (free acid)	[pending]	993.6 (free acid)	λ <sub>max</sub> 259 nm / ε 40500 / pH 7	T 064

**Name:** 2'-, 2'', 2'''- Trideoxy- 2'-, 2'', 2'''- trifluoro- cyclic triadenosine monophosphate, sodium salt  
Syn.: cFA3

**Description:** 2',2'',2'''-Tri-F-c-tridAMP is a cyclic nucleotide in which three 2'-deoxy-2'-fluoro-modified 5'-AMP units are interconnected via 3'-5' phosphodiester bonds to form a cyclic structure.

**Properties:** Cyclic oligoadenylates such as c-triAMP (Biolog Cat. No. C 362) were found to be novel bacterial second messengers involved in the Type III CRISPR-Cas-associated detection and degradation of invasive genetic elements in many prokaryotes. 2',2'',2'''-Tri-F-c-tridAMP is an analogue of c-triAMP containing fluoro modifications in all ribose 2'-positions which are expected to make the compound resistant to degradation by nucleases.

**Specification:** Crystallized or lyophilized sodium salt. Please keep in mind that equal concentrations of the compound may look different in volume due to sensitivity of the lyophilized form to humidity. The compound can even contract to small volume droplets. Normally the product is located in the conical bottom of the tube. Micromolar quantities are determined by UV at λ<sub>max</sub>.

**Purity:** Typical analysis is better than 95% (HPLC / UV / 259 nm). The product is not sterile and has not been tested for endotoxins.

**Solubility:** 2',2'',2'''-Tri-F-c-tridAMP is soluble in water and aqueous buffers (≥ 2.7 mM, limits have not been determined). Please rinse tube walls carefully and preferably use ultrasonic or vortex to achieve total and uniform mixing. When opening the tube please make sure that no substance is lost within the cap.

**Stability and Storage:** 2',2'',2'''-Tri-F-c-tridAMP has sufficient stability at room temperature and does not need special care during handling or shipment. Nevertheless, we recommend that the compound should be stored in the freezer, for longer storage periods preferably in freeze-dried form.

**Toxicity and Safety:** Please keep in mind, that the *in vivo* properties of this compound are not sufficiently characterized up to now. Avoid contact with eyes and skin or ingestion and allow only trained personnel to handle the product. Our products are designed, developed and sold for research purposes only! They are intended for *in vitro* and nonhuman *in vivo* laboratory applications. Any other use requires approval of health authorities.

**Not for drug, household or related uses!**

**Selected References for 2',2'',2'''-Tri-F-c-tridAMP:** 2',2'',2'''-Tri-F-c-tridAMP is a new product and there are currently no references available.

**Selected References for the Parent Compound c-triAMP:**

Lau, R.K.; Ye, Q.; Birkholz, E.A.; Berg, K.R.; Patel, L.; Mathews, I.T.; Watrous, J.D.; Ego, K.; Whiteley, A.T.; Lowey, B.; Mekalanos, J.J.; Kranzusch, P.J.; Jain, M.; Pogliano, J.; Corbett, K.D., *Mol. Cell*, **77**, 723 - 733 (2020): "Structure and Mechanism of a Cyclic Trinucleotide-Activated Bacterial Endonuclease Mediating Bacteriophage Immunity"

Ye, Q.; Lau, R.K.; Mathews, I.T.; Birkholz, E.A.; Watrous, J.D.; Azimi, C.S.; Pogliano, J.; Jain, M.; Corbett, K.D., *Mol. Cell*, **77**, 709 - 722 (2020): "HORMA Domain Proteins and a Trip13-like ATPase Regulate Bacterial cGAS-like Enzymes to Mediate Bacteriophage Immunity"

Smalakyte, D.; Kazlauskienė, M.; Havelund, J.F.; Rukšėnaitė, A.; Rimaite, A.; Tamulaitienė, G.; Færgeman, N.J.; Tamulaitis, G.; Siksnys, V., *Nucleic Acids Res.*, **48**, 9204 - 9217 (2020): "Type III-A CRISPR-associated Protein Csm6 Degrades Cyclic Hexa-adenylate Activator Using Both CARF and HEPN Domains"

Kazlauskienė, M.; Kostiuk, G.; Venclovas, Č.; Tamulaitis, G.; Siksnys, V., *Science*, **357**, 605 - 609 (2017): "A Cyclic Oligonucleotide Signaling Pathway in Type III CRISPR-Cas Systems"